YS600-YS800

**MANUAL DE USUARIO
USER MANUAL
BEDIENUNGSANLEITUNG**

Cavitación

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USER MANUAL YS600-YS800

English

Congratulations, first of all we would like to thank you for purchasing CELL CONTOUR and for the trust you have placed in us. By means of this handy User Manual and following its simple instructions you will achieve results comparable to those obtained from the hands of a professional beautician.

We strongly recommend that you read this User Manual thoroughly, reading and understanding it will help you use CELL CONTOUR effectively and safely

CONTENTS

Chapter 1

- Contents of the box
- Assembling CELL CONTOUR

Chapter 2

- Basic characteristics of CELL CONTOUR
- Operating instructions

Chapter 3

- Application modes

Chapter 4

- Correct, safe use of CELL CONTOUR

Chapter 5

- Frequently asked questions about CELL CONTOUR

Chapter 6

- Important safety instructions
- Warranty
- Certificates
- Parts list
- Technical specifications


Chapter 1
CONTENTS OF THE BOX

Inside this box you will find everything you need to enjoy CELL CONTOUR from now onward. Carefully check that the pack contains all of the elements on the following list.

YS600

- 10. Main unit
- 11. Hand unit support with screws
- 12. Hand unit 1
- 13. Mains cable
- 14. Cavitation gel
- 15. Toning serum
- 16. Anti-cellulite serum
- 17. Post-cavitation gel
- 18. User manual


Fig.1

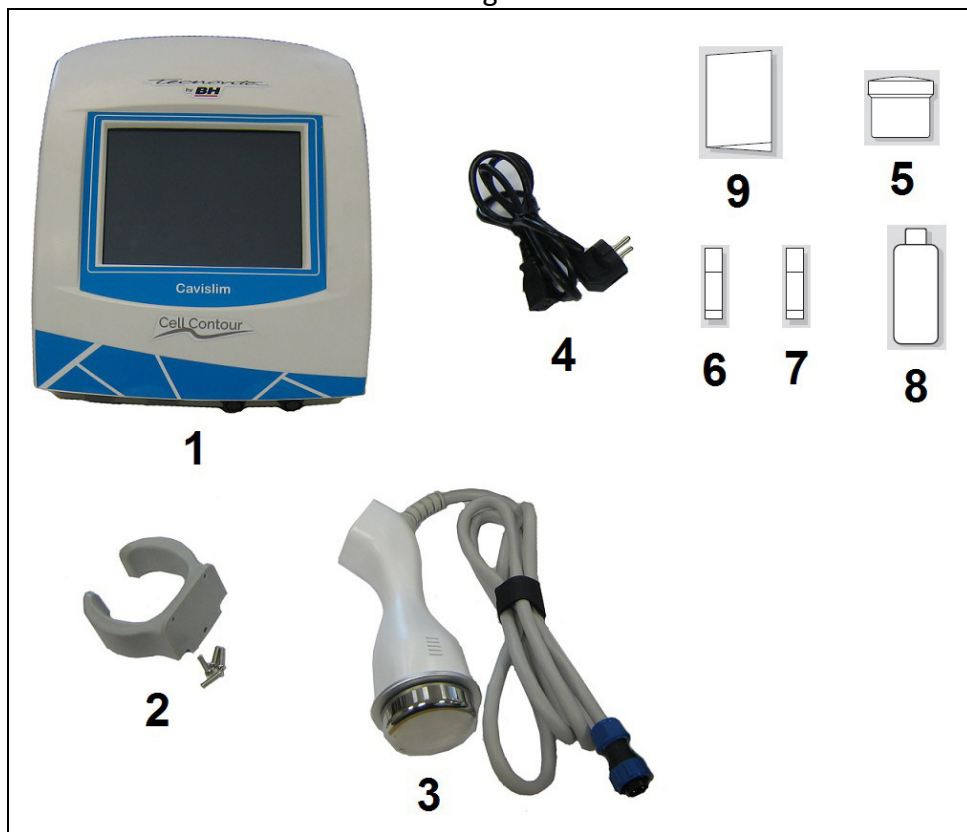YS800

1. Main unit
2. Hand unit support with screws
3. Hand unit 1
4. Hand unit 2
5. Mains cable
6. Cavitation gel
7. Toning serum
8. Anti-cellulite serum
9. Post-cavitation gel
10. User manual

Fig.2

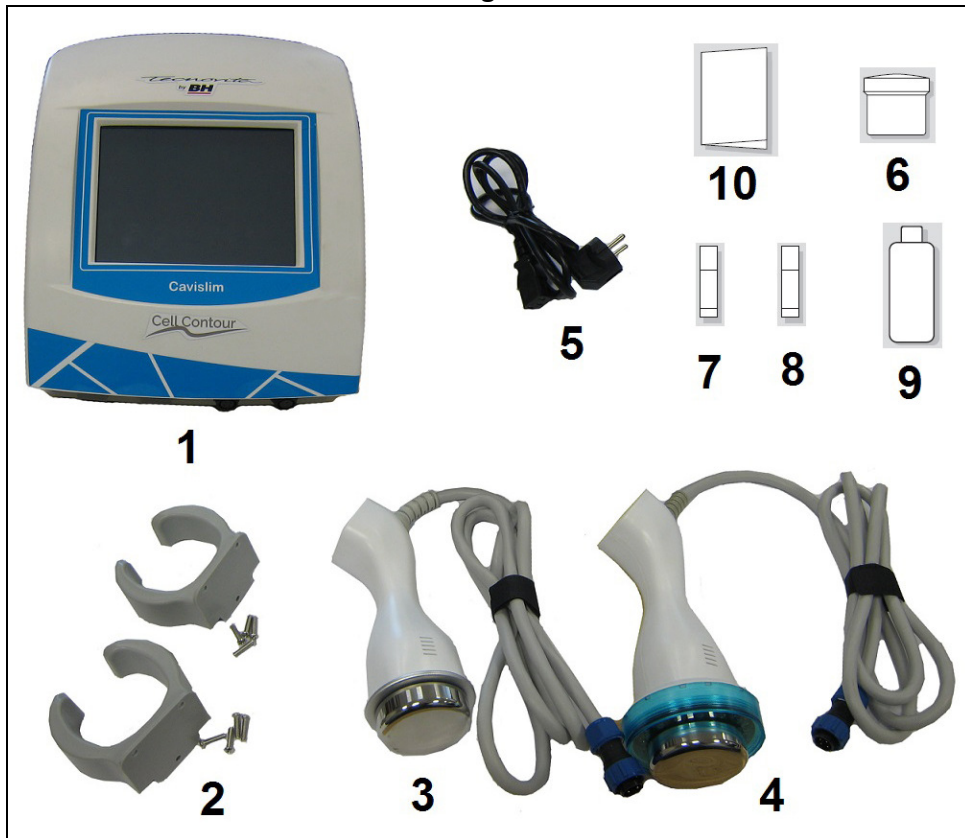
ASSEMBLING CELL CONTOUR

First we are going to fit the handset support onto the main unit. To do this, place the main unit on a table and position the handset support by inserting the screws as shown in Fig. 3.

Using the fixed wrench tighten the screws securely so that the two parts fit together perfectly.

Fig.3

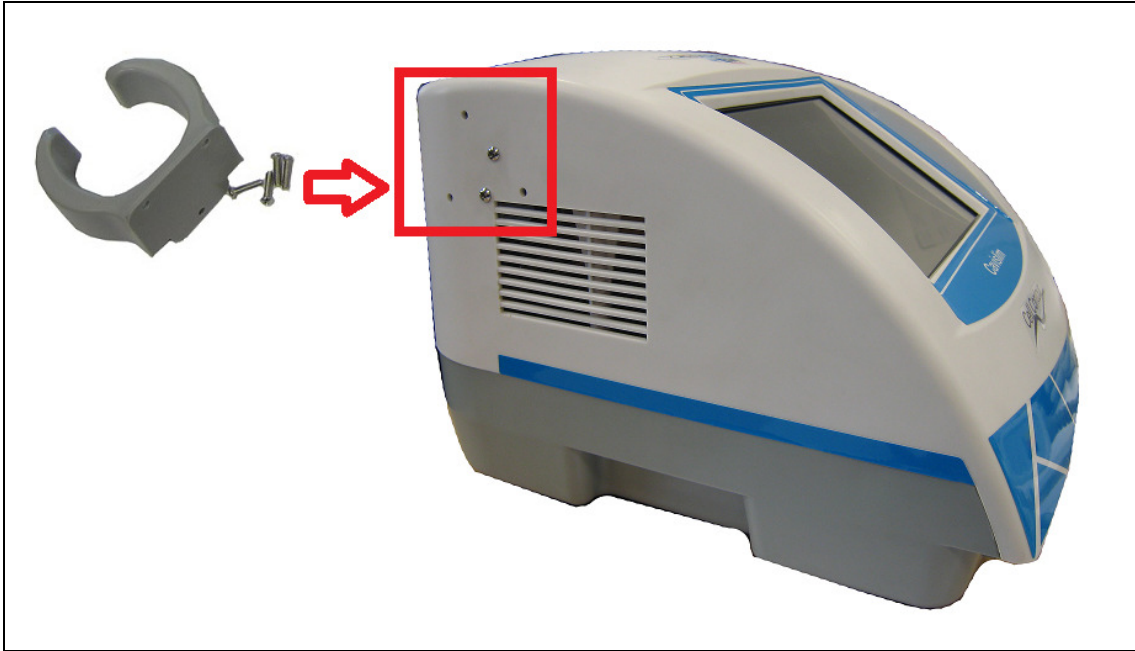Chapter 2

BASIC CHARACTERISTICS OF CELL CONTOUR

CELL CONTOUR is a hi-tech device that allows us to selectively and efficiently remove the build-up of fat (cellulite, orange skin syndrome) from our body.

The technology used is that of ultrasound emission which produces the physical phenomenon known as cavitation, causing the removal of fat cells.

This technology with proven results is widely used in beauty salons all over the world in order to fight cellulite painlessly and non-invasively.

The intensity function controls how deep the sound waves penetrate your skin. The greater the intensity, the deeper the waves reach.

The time function allows you to choose how long you wish to work with CELL CONTOUR. For reasons of safety the time function allows you to choose up to 20 minutes which corresponds to the maximum recommended time for working a specific area.

CELL CONTOUR also lets you choose the work mode, helping you select the time interval during which the waves reach your skin.

In the auto-programs you can adjust the intensity, the rest of the features is already set up.

Characteristics for model CELL CONTOUR Ref. YS600

Time: Adjustable from 1 to 20 minutes.
Intensity: Adjustable 7 intensity levels (3,5-6W/cm²)
Mode: Pulsed or continuous
Programs: 4 Pre-set programs
 1 Manual program
Handset: 1

Characteristics for model CELL CONTOUR Ref. YS800

Time: Adjustable from 1 to 20 minutes.
Intensity: Adjustable 7 intensity levels (3,5-6W/cm²)
Mode: Pulsed or continuous
Programs: 4 Pre-set programs
 1 Manual program
Handset: 2

The YS800 model has an additional handset that combines the technology of ultrasound with infrared technology.

How does the infrared light affect our skin?

The infrared light helps promote the production of collagen and elastin, making the skin firmer and more flexible. It promotes blood circulation, which generates new cells and makes the skin look healthier and rejuvenated.

CELL CONTOUR OPERATING INSTRUCTIONS


Precautions: Consult your doctor before beginning to use this appliance. This warning is especially important for people with health problems. Use by those under 18 years of age is not recommended. Keep this appliance out of children's reach.

Starting

Connect the hand unit (two hand units for YS800) to the front of the machine.


Connect the hand unit by matching the numbers of the male pins with the numbers of the female pins (Fig. 4).

Fig.4

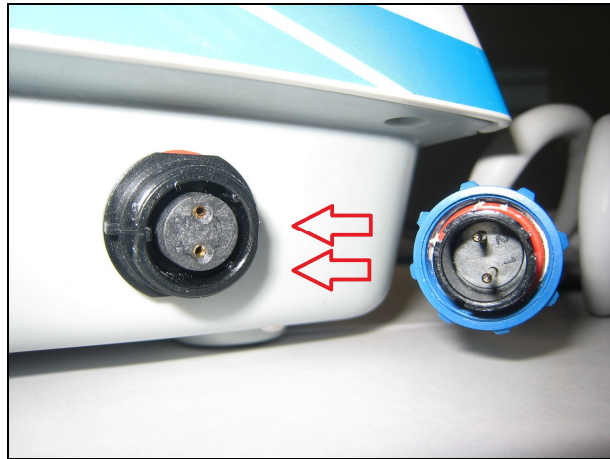Then secure the connection by turning the plastic thread in the clockwise direction (Fig.5).

Fig.5

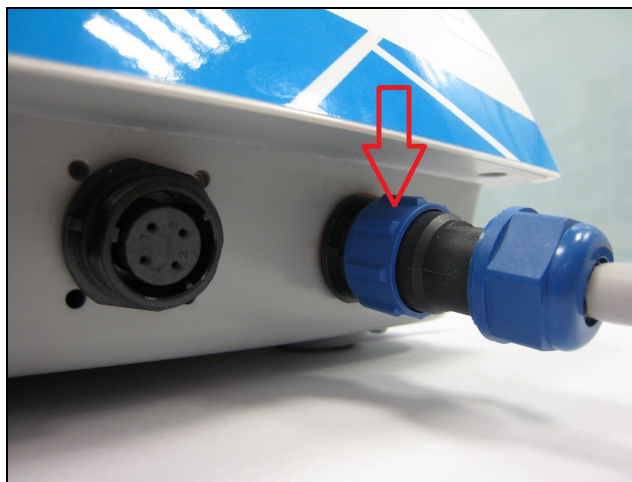Check that all the connections are made correctly before connecting it to the mains supply.


Finally, plug one end of the mains cable into the socket on the machine (Fig.6) and the other end into the electric wall socket (220V).

Fig.6

Switch the On/Off switch to position (I) to start the machine (Fig.7).

Fig.7

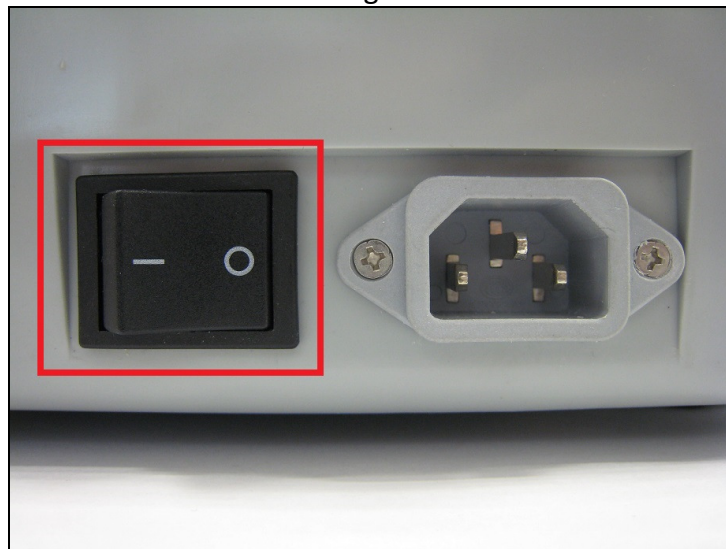
Monitor elements

CELL CONTOUR has a latest technology touch screen.


YS800

Fig.8

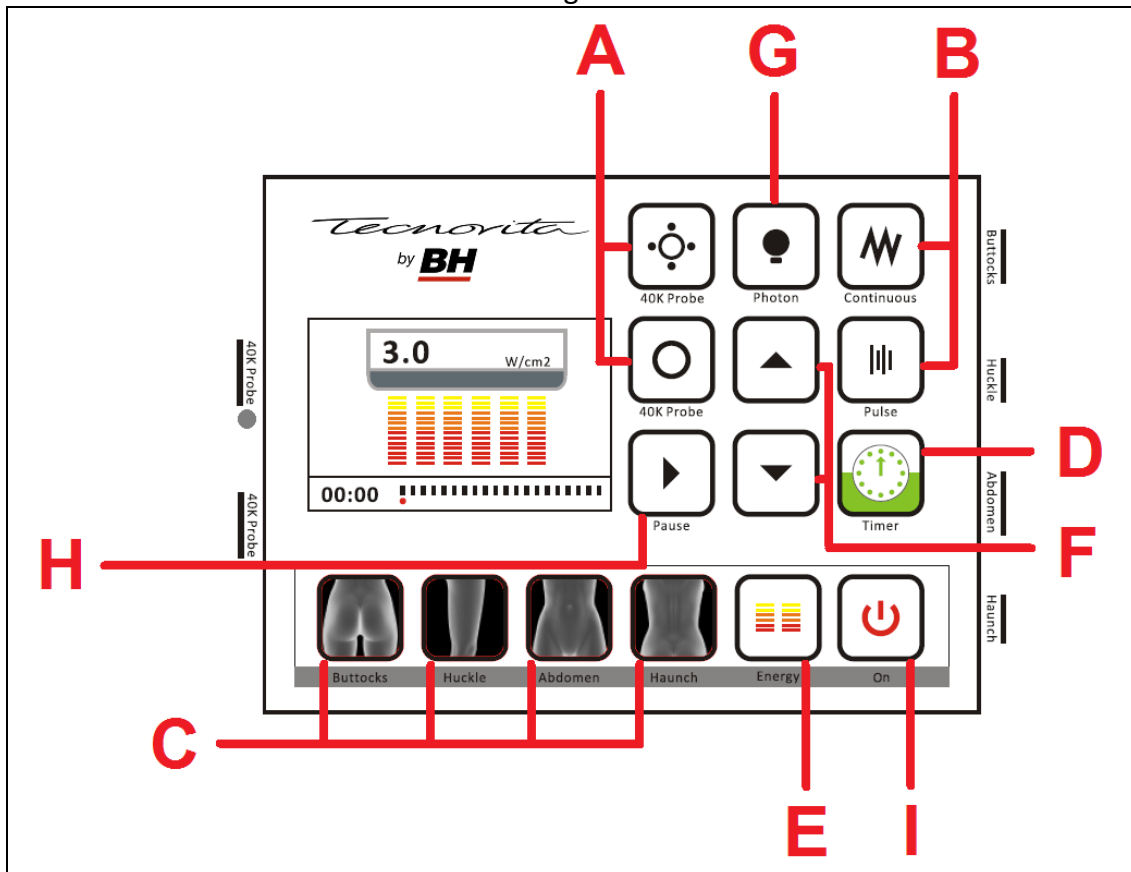- A – Selection of hand unit in manual program
- B – Selection of the mode (continuous or pulsed)
- C – Selection of the pre-set programs (buttocks, huckle, abdomen, haunch)
- D – Time adjustment
- E – Intensity adjustment
- F.- Time/intensity adjustment after selecting time/intensity
- G – Turn on/turn off of the infrared light function
- H – Pause during a session
- I – Start/stop of a session

Fig.9

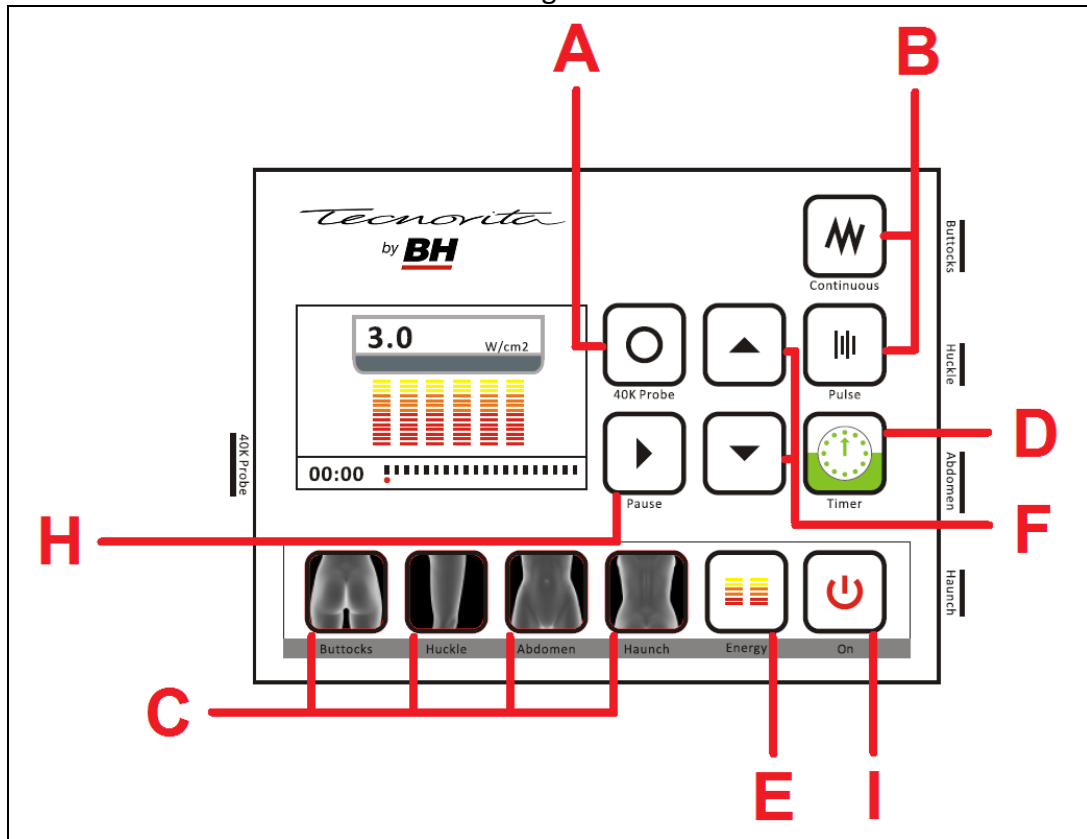- A – Selection of hand unit in manual program
- B – Selection of the mode (continuous or pulsed)
- C – Selection of the pre-set programs (buttocks, huckle, abdomen, haunch)
- D – Time adjustment
- E – Intensity adjustment
- F.- Time/intensity adjustment after selecting time/intensity
- H – Pause during a session
- I – Start/stop of a session

* If you are inside the time/intensity adjustment in the manual program or intensity adjustment in pre-set program, you have to click on the time/intensity icon to exit the adjustment first and then you can select other functions.

Cavitation session

There are 3 settings for the massage for YS600 and 4 settings for YS800:

- **Time:** measures the duration of the session in minutes and seconds. You can select between 1 to 20 minutes (in steps of 1 minute).

- **Work mode:** Pulse type (impulse emission frequency).

CELL CONTOUR allows you to choose between pulsed and continuous mode.

- **Intensity:** establishes the intensity level for the session.

CELL CONTOUR allows you to choose different intensity levels (7 levels).

- **Handset (for YS800):**

CELL CONTOUR ADVANCE allows you to choose between the handset with ultrasound function and the handset with ultrasound and infrared function to improve the results of the sessions.

Beginning a session

Once you have plugged the machine into the mains supply and switched it to position (I), the monitor will show a menu to select the language. Click on the flag of the desired language. You will enter the screen of the selection of programs (Fig.8).

Fig.10 - YS800


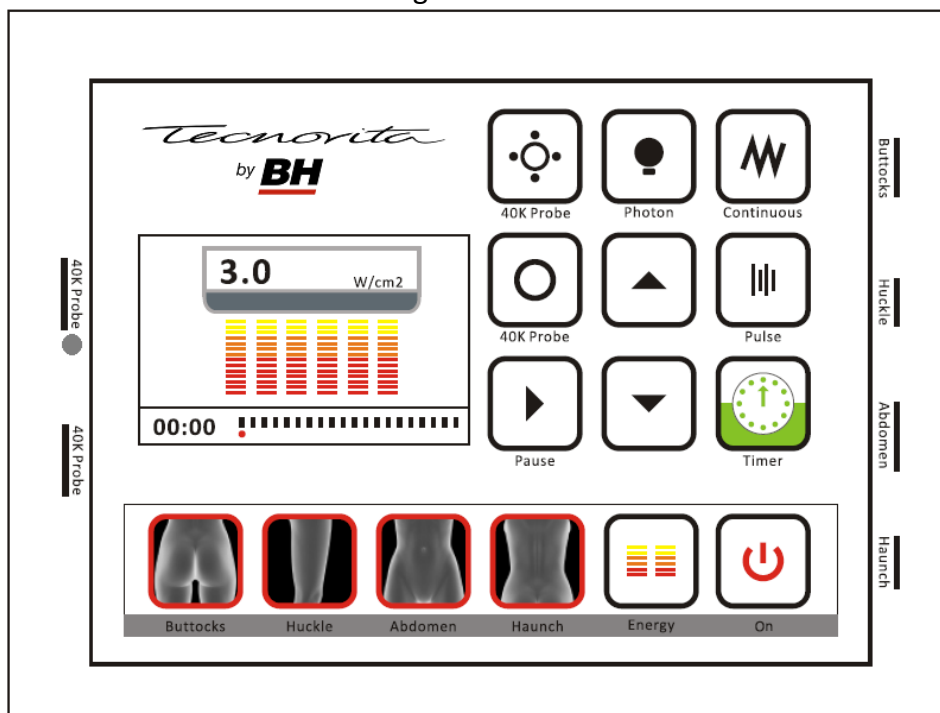
Fig.11 – YS600

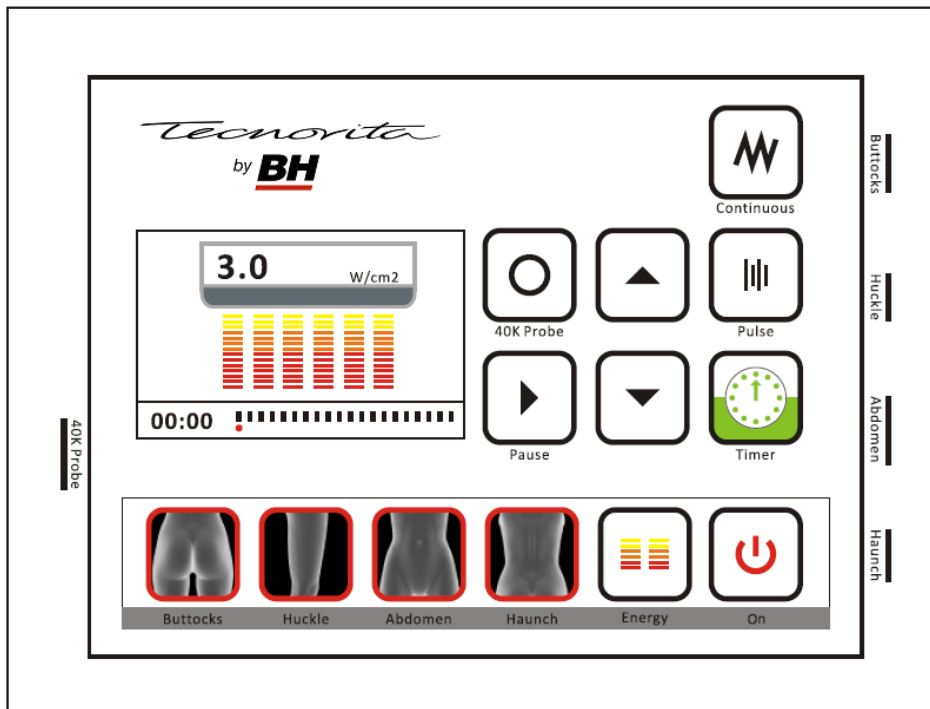Always apply an abundant amount of gel to the area that you wish to work. If you need to purchase gel or any of our creams you can find them in the shop where you purchased CELL CONTOUR or at our website: www.tecnovita.eu

Pre-set profiles

CELL CONTOUR has 4 pre-set profiles designed specifically for these body areas: buttocks, huckle, abdomen, haunch. Click on the icon of the area you want to work. You will see that the pre-set values of the program will be marked and the hand unit will be selected. You can modify the intensity by pressing the intensity icon (E) and then the adjustment buttons (F). Once the body part has been selected and the intensity adjusted (in case you want to modify it), click on the icon for start / end of the session (I) to start the cavitation session. Begin to treat the area by sliding the hand unit slowly and without pressure, as shown in the pictures in the chapter of Application Modes. You can modify the intensity during the cavitation session by pressing on intensity adjustment icon (E).

The unit will beep every 30 seconds. If you want to pause the session, click on the pause icon (H). To return to the session, press the pause icon (H) again.

When the time counter reaches **0.00** the unit will stop and emit a series of beeps to tell you that the selected session has finished.

Characteristics of the pre-set profiles:

Body part	Infrared (YS800)	LED	Power/Intensity (W/cm2)	Time (min)	Mode
Buttocks	Sí		4,5	10	Continuo
Huckle	No		5	12	Continuo
Abdomen	Sí		5,5	8	Continuo
Haunch	No		4	10	Continuo

Manual program

Select the hand unit by clicking on the hand unit icon (A). To turn on / turn off the infrared LED function (for the YS800), click on the icon (G) once you have selected the infrared hand unit. Select the work mode by clicking on (B). Adjust the time by clicking on the time icon (D) and then on the adjustment buttons (F). Adjust the intensity by pressing on the intensity icon (E) and then on the adjustment buttons (F). Click on the start / stop icon (I) to start the session. Begin to treat the area by sliding the hand unit slowly and without pressure, as shown in the pictures in the chapter of Application modes. You can modify the intensity during the cavitation session by pressing on intensity adjustment (E).

The unit will beep every 30 seconds. If you want to pause the session, click on the pause icon (H). To return to the session, press the pause icon (H) again.

When the time counter reaches **0.00** the unit will stop and emit a series of beeps to tell you that the selected session has finished.

Chapter 3

CELL CONTOUR APPLICATION MODES

Always apply an abundant amount of gel to the area that you wish to work and begin to treat it by sliding the handset slowly and without pressing, as shown in the pictures. If you need to purchase gel or any of our preparation creams you can find them in the shop where you purchased CELL CONTOUR or at our website: www.tecnovita.eu


In each of the following sections you will find a table with application protocols for the manual program. Before beginning to use the appliance, you must determine the level of fat in each area to be treated.

- A Large accumulation of high density fat
- B Large accumulation of medium density fat
- C Medium accumulation of medium density fat
- D Small accumulation of fat
- E orange peel syndrome
- F Maintenance session

CENTRAL ABDOMEN

Remember, never apply CELL CONTOUR to the navel, breasts, the thorax area or the genitals.

Apply an abundant amount of gel on your abdomen and begin to massage the area by sliding the handset slowly and without pressing, as shown in the diagrams.

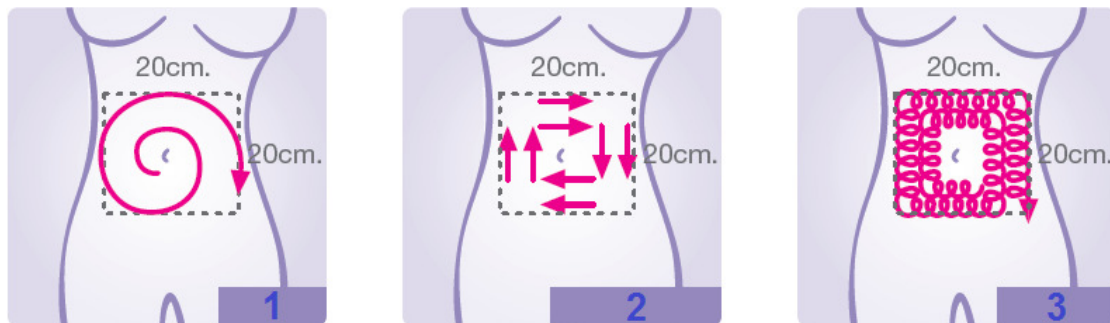For best results we recommend applying CELL CONTOUR by combining different massage techniques (1-spiral,2-crossways,3-circular) in different sessions.


We recommend lying down for a more comfortable CELL CONTOUR application.

Level	Mode	LED (YS800)	Intensity	Time
A	Continuous	No	5,5	10'
B	Continuous	No	5,5	10'
C	Continuous	Yes	5,5	8'
D	Pulsed	Yes	4	8'

E	Pulsed	Yes	4	8'
F	Pulsed	Yes	3,5	8'

SIDES

Apply gel over the area to be worked using CELL CONTOUR

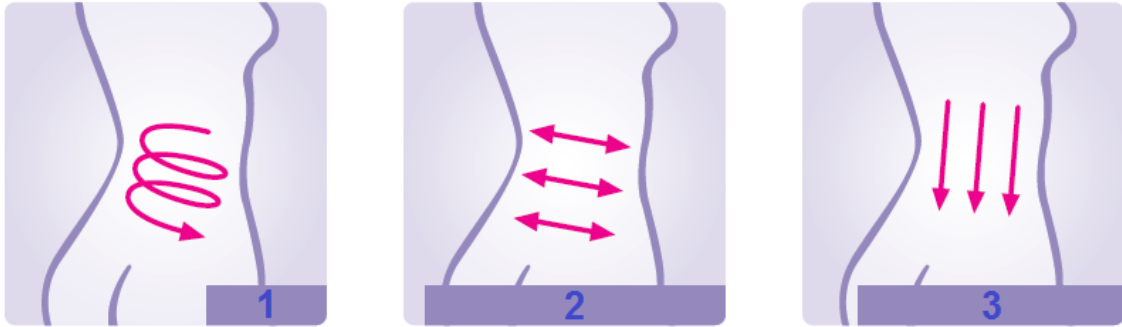
Circular (1) and longitudinal (2,3) massage movements should be applied from top to bottom as shown in the diagrams. It is important to massage in the correct direction, as this will improve drainage.

Combine the various massage techniques during different sessions.

Level	Mode	LED (YS800)	Intensity	Time
A	Continuous	No	5	10'
B	Continuous	No	5	10'
C	Pulsed	Yes	4,5	7'
D	Pulsed	Yes	4	7'
E	Pulsed	Yes	4	5'
F	Pulsed	Yes	3,5	5'

HIPS

Apply gel over the area to be worked using CELL CONTOUR.

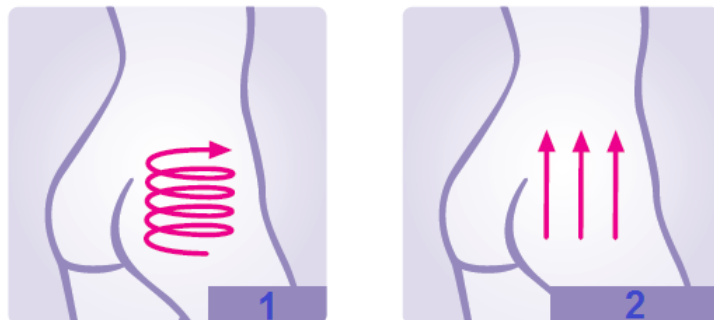
Circular (1) and longitudinal (2) massage movements should be applied from the bottom upward as shown in the diagrams.

Level	Mode	LED (YS800)	Intensity	Time
A	Continuous	No	5	12'
B	Continuous	No	4,5	10'
C	Continuous	No	4	7'
D	Pulsed	Yes	3,5	5'
E	Pulsed	Yes	3,5	5'
F	Pulsed	Yes	3	5'

BUTTOCKS

Apply gel over the area to be worked using CELL CONTOUR.

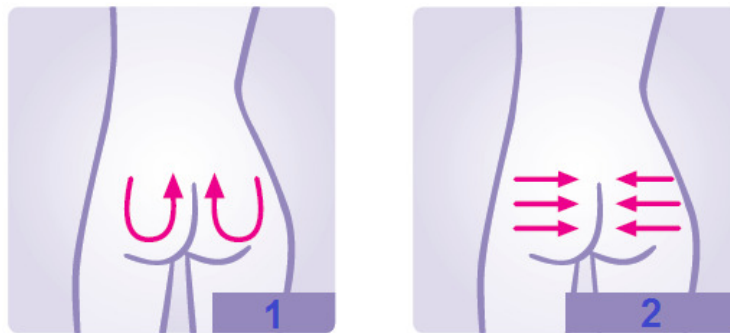
Both longitudinal (2) and circular (1) massage should be done in the direction of the arrows shown in the diagrams.

Level	Mode	LED (YS800)	Intensity	Time
A	Continuous	No	5,5	12'
B	Continuous	No	5	12'
C	Continuous	Yes	4,5	10'
D	Pulsed	Yes	2	10'
E	Pulsed	Yes	2	10'
F	Pulsed	Yes	2	8'

INNER THIGH

Remember, never apply CELL CONTOUR from the knee to the foot, there is no significant accumulation of fat in this area.

Apply gel over the area to be worked using CELL CONTOUR.

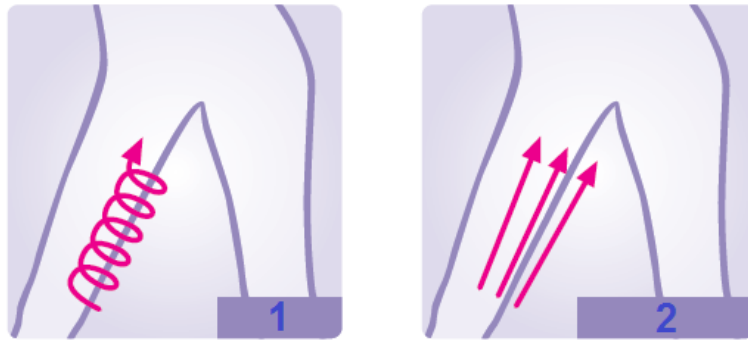
Both the longitudinal (2) and circular (1) massage should be done from the bottom upward as shown in the diagrams. Follow these simple instructions to help you get the best results.

Level	Mode	LED (YS800)	Intensity	Time
A	Pulsed	No	4	10'
B	Pulsed	No	4	10'
C	Pulsed	Yes	3,5	7'
D	Pulsed	Yes	3	7'
E	Pulsed	Yes	3	7'
F	Pulsed	Yes	3	7'

INNER KNEE

Remember not to apply CAVISILIM on the knee joints.

Apply gel over the area to be worked using CELL CONTOUR.

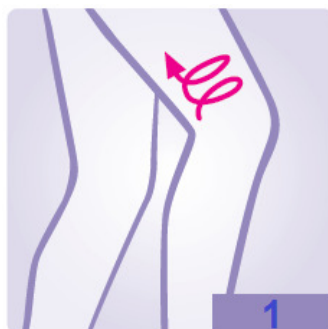
Apply a light circular (1) massage to this area, as shown in the diagram.

Level	Mode	LED (YS800)	Intensity	Time
A	Pulsed	No	3,5	10'
B	Pulsed	No	3,5	10'
C	Pulsed	Yes	3	7'
D	Pulsed	Yes	3	7'
E	Pulsed	Yes	3	7'
F	Pulsed	Yes	3	7'

ARMS

Remember, never apply CELL CONTOUR on the hands, forearms or joints.

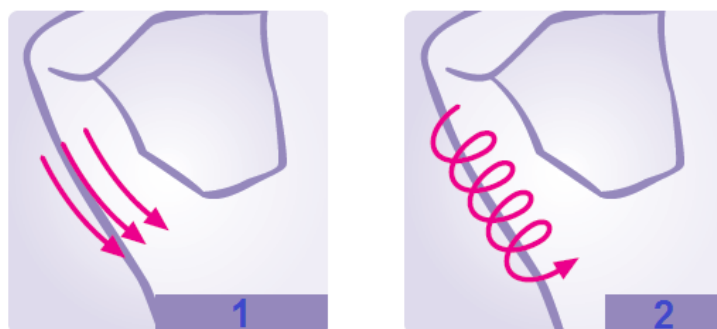Put a generous amount of gel on the area to be worked and this will improve the effectiveness of the treatment and prevent damaging the machine.

Perform the massage, sliding the handset slowly and without pressing. Massaging this area should be applied toward the body, as shown in the diagram.

Level	Mode	LED (YS800)	Intensity	Time
A	Pulsed	Yes	3,5	7'
B	Pulsed	Yes	3,5	5'
C	Pulsed	Yes	3	5'
D	Pulsed	Yes	3	3'
E	Pulsed	Yes	3	3'
F	Pulsed	Yes	3	5'

Chapter 4

Correct, safe use of CELL CONTOUR

For correct and safe use of CELL CONTOUR, you should keep the following application times in mind:

- Per work area: 20 minutes maximum.
By "work area" we mean each one of the body parts where CELL CONTOUR can be applied. The various work areas are: arms, abdomen, waist, hips, inner thighs and inner knees.
- Per session: 30 minutes maximum.
By session we mean using CELL CONTOUR on two different work areas (15 minutes per area) in the same day.
E.g.: In a single session you can apply CELL CONTOUR for 15 minutes on your sides and 15 minutes on your arms. But, for example, you must NEVER perform a 30 minute application session just on your sides in the same day.
- Between sessions: 72 hours minimum.
By "between sessions" we mean the rest time between each session. *E.g.: In one week the session schedule (maximum) would be: Monday, Thursday and Sunday.*
- Every 20 sessions: You should rest for 3 weeks. After this period you may decide to repeat the first treatment or start the maintenance program.
- Maintenance sessions: This depends on the individual, but one session per month should suffice.

VERY IMPORTANT

***** Never exceed the application times (20 minutes per area, 30 minutes total in one day).** Your body does not know "where" you are eliminating the fat from, so a session must never include 20 minutes on your sides, 20 minutes on your legs, 20 minutes on your arms, etc. REMEMBER: the maximum time during one day, adding all the parts of the body being treated, must be no more than 30 minutes.

***** And, above all, respect the rest period of 72 hours between sessions.** This is the time that your body needs to eliminate all the fat dissolved during the session through the lymphatic system.

Chapter 5

Frequently asked questions

What is cavitation?

Cavitation is a non-invasive (surgery free) body shaping technique that reduces volume by eliminating cellulite and localised fat deposits in the abdomen, hips, flanks, arms, inner thighs and knees.

What are the advantages over other methods?

The main advantage is that we can choose the specific areas of the body that we want to work. It could be said that it is a technology which allows us to concentrate on cellulite selectively.

Is cavitation painful?

Use of CELL CONTOUR should never cause pain.

Is it normal to hear a slight buzzing in your ears while using CELL CONTOUR?

Yes, some users may detect a slight buzzing sound during its application.

Can it be used on all parts of the body?

It cannot be applied to all parts of the body. Under no circumstances can it be applied to: the head or face, neck, breasts, the thorax, navel, forearm, joints, genitals, back, ankles or lymph glands.

Does CELL CONTOUR reduce weight and volume?

CELL CONTOUR eliminates clusters of fat, which have weight; consequently your body loses volume.

Is CELL CONTOUR effective on orange peel syndrome?

Yes, CELL CONTOUR is effective on orange peel syndrome because it destroys adipose cells that cause fat storage.

Can CELL CONTOUR be used to remove wrinkles?

No, CELL CONTOUR has not been designed to remove wrinkles.

Can CELL CONTOUR be used on the face?

CELL CONTOUR must NEVER be used on the face, eyes or neck.

After how many sessions will I begin to see results?

You should notice results after the first five sessions, however results may vary depending on the subject.

How long should I wait between sessions?

It is advisable to wait 72 hours between sessions.

Does cellulite reappear after using CELL CONTOUR?

Although CELL CONTOUR ensures effective and long-lasting results, we cannot guarantee against cellulite reappearing. To prevent the appearance of new cellulite we recommend a balanced diet, regular physical exercise and maintenance sessions using CELL CONTOUR.

Are there any side effects?

If you are healthy, CELL CONTOUR has no reported side effects; nevertheless in some people it may produce a reddening of the skin that can last for about 2 hours.

Is it necessary to combine CELL CONTOUR with a special diet?

No, it is not necessary to follow any special diet, but a balanced diet and proper hydration is highly advisable.

How much water should I drink?

It is essential to begin hydrating yourself two hours before applying CELL CONTOUR. It is also advisable to drink more water after finishing a session. The recommended quantity is at least 2 litres during these periods, as it helps to eliminate toxins.

Is it necessary to combine CELL CONTOUR with a sports activity?

No, this is not strictly necessary, but at least 30 minutes of physical activity will enhance the results as it helps to drain fluids.

Which sports activities can supplement CELL CONTOUR?

The vibration platform is the ideal complement due to its toning and draining qualities. But if you do not have access to a platform, then walking, exercise bike, swimming or any other aerobic exercise for at least 30 minutes is very helpful.

How can I enhance the effects of CELL CONTOUR?

By combining a balanced diet with regular physical exercise and drinking plenty of water. Added to which, drainage massage can also be very useful.

What are the anti-cellulite and toning gels for? When do I use them?

1. Anti-cellulite gel: Contains active ingredients that help eliminate cellulite.
2. Toning gel: Contains active ingredients that make the treated area firmer.

The application sequence is as follows: during the first few days of treatment use the anti-cellulite gel and once you begin to notice results, use the toning gel This is applied before the Cavitation gel, as with any other serum. The area should be massaged until it is well absorbed and then apply plenty of cavitation gel to start the session

Can I use CELL CONTOUR without applying gel?

NEVER, this is because the gel acts as a conductor for the ultrasonic sound waves with your body. If you do not use gel then the CELL CONTOUR application loses its effectiveness and the machine can be severely damaged.

Can I use any type of gel?

We recommend using only TECNOVITA by BH gel, as the use of any other type may damage the machine.

When is it inadvisable to use CELL CONTOUR?

- If you are pregnant or think you might be.
- If you are breast-feeding.
- Undergoing fertility treatments.
- If you have any kind of heart disease.
- If you suffer from skin problems.
- Osteoporosis.
- Fibromyalgia.
- Malignant tumours.
- Hypertension.
- High cholesterol levels.
- Kidney problems
- Physical shortcomings
- Metal implants or prostheses.
- Physical malformations.
- Epilepsy.
- Long-term treatment with medications: anti-inflammatories, anticoagulants, antibiotics or antihistamines.
- Depression.
- Psychiatric or psychological treatment.
- Liver problems (fatty liver)
- Irritable bowel syndrome.
- Diabetes

Can I use CELL CONTOUR if I am pregnant?

You must never use CELL CONTOUR if you are pregnant or think you might be.

Can CELL CONTOUR be used during lactation?

Do not use it during lactation.

Can CELL CONTOUR be used on the breasts?

NEVER apply it to the breasts.

Can I use CELL CONTOUR during menstruation?

Use on the abdominal area is not advisable during menstruation as it may cause increased bleeding. However, it can be applied in other areas.

Can I use CELL CONTOUR if I have an IUD?

Yes.

Can CELL CONTOUR affect my cholesterol levels?

Use of CELL CONTOUR may increase your cholesterol levels.

Can CELL CONTOUR be used with skin conditions?

It should not be used if you have any skin condition. Consult you doctor.

Can I use CELL CONTOUR if I have ear implants?

Do not use it if you have implants fitted.

Can I use CELL CONTOUR when I am bathing or showering?

Under no circumstances should the handset or the machine come into contact with water due to the risk of electrocution.

How do I clean CELL CONTOUR and the handset?

Never immerse CELL CONTOUR or the handset in water. Use a cloth slightly dampened with water.

Chapter 6

IMPORTANT SAFETY ADVICE

This product is supplied with a mains cable that has an earth wire and a plug with an earth pin.

In the event that the appliance malfunctions then the earth connection will redirect the electrical current away from the machine and reduce the risk of an electrical shock.

WARNING!: If the equipment is not earthed correctly then there is a risk of electrical shock.

Do not modify the plug supplied with the equipment, if it does not fit the wall socket then get an electrician to install an electrical socket that does.

This product should be used on a 220 volt circuit.

Follow these recommendations:

6. Never use accessories that are not recommended by the manufacturer.
7. Do not use the machine if the mains cable becomes damaged or worn.
8. Keep the mains cable away from hot surfaces.
9. Do not use the machine outdoors.
10. To unplug the machine, first put the switch to the «0» position and then unplug it from the wall socket.

Warranty

We suggest that you read through this document carefully in order to fully understand the cover that the manufacturer offers as well as your responsibilities as owner.

Guarantees

GENERAL INFORMATION

EXERCYCLE, S.L. would like to thank you for purchasing this product and the trust you have placed in us. Our products have been subjected to rigorous quality checks. Nevertheless, should you experience any problem with this appliance during the guarantee period then please read through these warranty conditions and then get in touch with your local authorised dealer.

This guarantee certificate has been drafted in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding on certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product under the following terms and conditions:

2. GUARANTEE PERIOD

The product covered by this guarantee certificate is protected for a period of TWO (2) YEARS as from the date of receiving the product, which will be taken as the date that appears on the sales receipt for the product or on the corresponding delivery docket should this be later. It is essential that this be presented in the event of possible repairs. The policy regarding the warranty period may vary between different countries. The stated warranty period is only applicable to Europe.

2. RIGHTS UNDER THE GUARANTEE

- Repair or replacement of the product which, in such instance, will be carried out within a reasonable period and without significant inconvenience to the consumer, taking into account of the nature of the goods and the purpose for which the consumer required them.
- Those cases in which the repair or replacement represents a disproportionate cost for the seller in comparison with other forms of remedy will be exempt assuming these do not prejudice the consumer.
- Repairs will be free for the consumer with regard to labour, materials and carriage.
- In the event that the machine is repaired or replaced with a new one, the guarantee period will be suspended temporarily until the new appliance has been delivered and will carry at least six months of guarantee, however under no circumstances will a new 2 year period be initiated.
- The buyer may, as is apparent, require (i) a reduction in price or proceed to a (ii) recession of the contract when the repair or replacement is not carried out within a reasonable time or with significant inconvenience to the buyer. The reduction in price (i) will be proportional to the difference between the value that the goods would have had if they had been delivered in accordance with the contract and the actual value that the goods were worth at the time of the said delivery. Recession of the contract (ii) will not proceed when the non-conformity is of minor importance.

3. THE GUARANTEE DOES NOT COVER

- Defects attributable to use other than that established in the User Manual.
- Professional use or commercial leasing or hiring of the product.
- Poor installation, poor assembly or breakage of components due to improper installation, (the User Manual gives detailed instructions on how to install and use the machine).
- Improper use of the product (use without applying gel, exceeding operating times, exposure to high moisture levels, etc.) or problems arising from poor maintenance (accumulation of dirt on any part of the product or its accessories, etc.) (check the User Manual for detailed instructions on how to

install and use the product).

- Repairs carried out by persons not authorised by EXERCYCLE, S.L. or appliances that have been opened or tampered with by persons not authorised by the manufacturer.
- Deviations from the standard product, such as light vibration, noise or breakage of components that do not affect the performance of the machine and which cannot be attributed to manufacturing faults.
- Scratches, scrapes, abrasions, knocks, falls, etc. caused by carelessness.
- Defects on the appliance caused by damage during transit.

4. CLAIMING AGAINST THE GUARANTEE

- Claims made for nonconformities will be limited to three years as from the date that the goods were delivered.
- The consumer must report the existence of the defect within a period of 2 months as from the date that the defect was discovered.
- This guarantee does not affect the consumer's statutory rights under any applicable law, unless it is incompatible with the application of other actions derived from compensation sought for hidden faults in the contract of sale.
- With regard to queries that may arise from the interpretation or application of this guarantee, these will be governed by the corresponding national law.

Any discrepancy arising from the interpretation, validity or application of this contract will be subject to the courts and tribunals of Vitoria, Spain.

Certificates

By virtue of this document, Exercycle S.L. declares that this product meets the essential requirements and other relevant provisions of Directives 2004/108/EC and 2006/95/EC

Parts list

To order replacement parts, please indicate:

The machine model | Corresponding part N° | Quantity

Example:

YS600	4	1
-------	---	---

YS600

5. Main unit
6. Handset support with screws
7. Handset 1
8. Mains cable

Fig.1

YS800

- 6. Main unit
- 7. Handset support with screws
- 8. Handset 1
- 9. Handset 2
- 10. Mains cable

Fig.2

Technical specifications

Net weight: 7,5 kg

Dimensions: 42 cm x 35 cm x 22 cm.

Ultrasound: 40 Khz

Power: 60W

Power supply: Voltage: AC 220V

Date of purchase

Purchased in

Buyer information

Dealer stamp

Please keep the receipt as proof of purchase

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